
Y-Connector


OKAY II

Y-Connector Set



Connecting
design,
function and
performance

Features

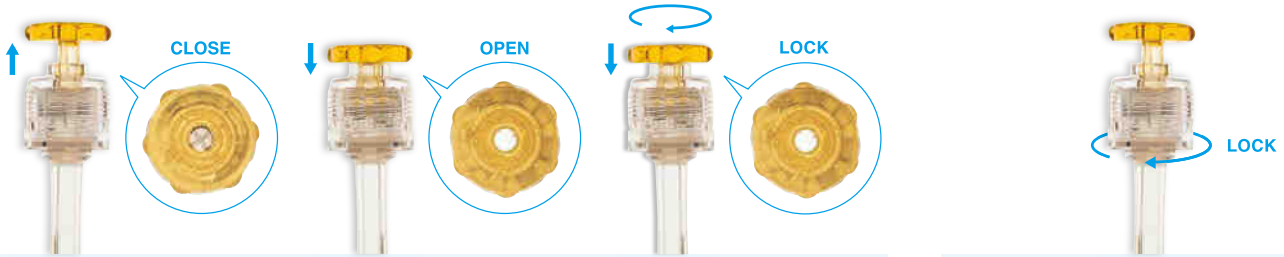

Hemostatic valve

The hemostatic valve can be opened by pushing the opener, and closed by releasing it.

Once the hemostatic valve has been opened, a clockwise quarter turn locks the valve in the open position.

Fixed valve

Valve to fix devices and prevent loss of pressure.

Large lumen

Nominal ID: 10F, 3.33mm

Allows for use in combination with various products for a variety of PCI strategies.

Superior design

Facilitates single-handed operation.


Ease-of-use

Hemostatic valve features a locking mechanism.

Minimizes blood loss

Two valves (hemostatic and fixed) greatly reduce blood loss during the procedure.

Accessories

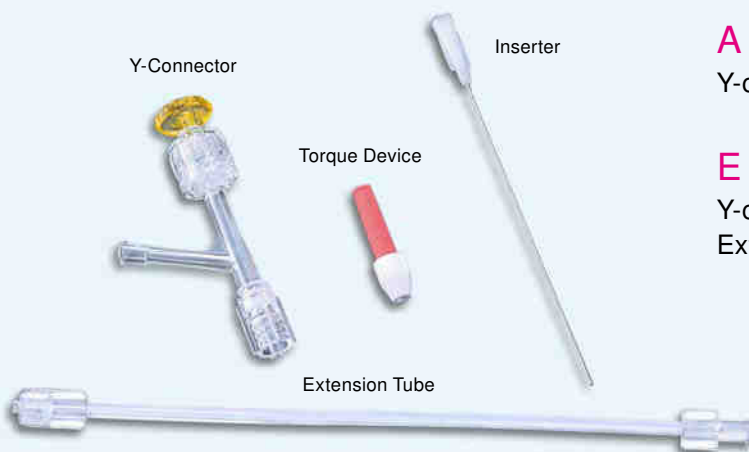

A set (YOK0A)

Y-connector, Inserter, Torque Device

E set (YOK0E)

Y-connector, Inserter, Torque Device, Extension tube

Manufacturer



<http://www.goodmankk.com/english/>

▶ Please consult your sales representative regarding stock availability.
▶ GOODTEC Y-Connector Set is a trademark of Goodman Co., Ltd.
▶ Please refer to the Instructions For Use for further details regarding products use and specifications.



OKAY II-1708-E6